


DELICATE PROCEDURES DEMAND PRECISION

The MedOne line of Auricle™ instruments deliver the ultimate in surgeon control for delicate procedures. Squeeze style mechanism and rotatable instrument tip enables precise control and comfortable access to the microsurgical site. Ideal for middle ear surgery, neurotology, and neurosurgery/skull base surgery.

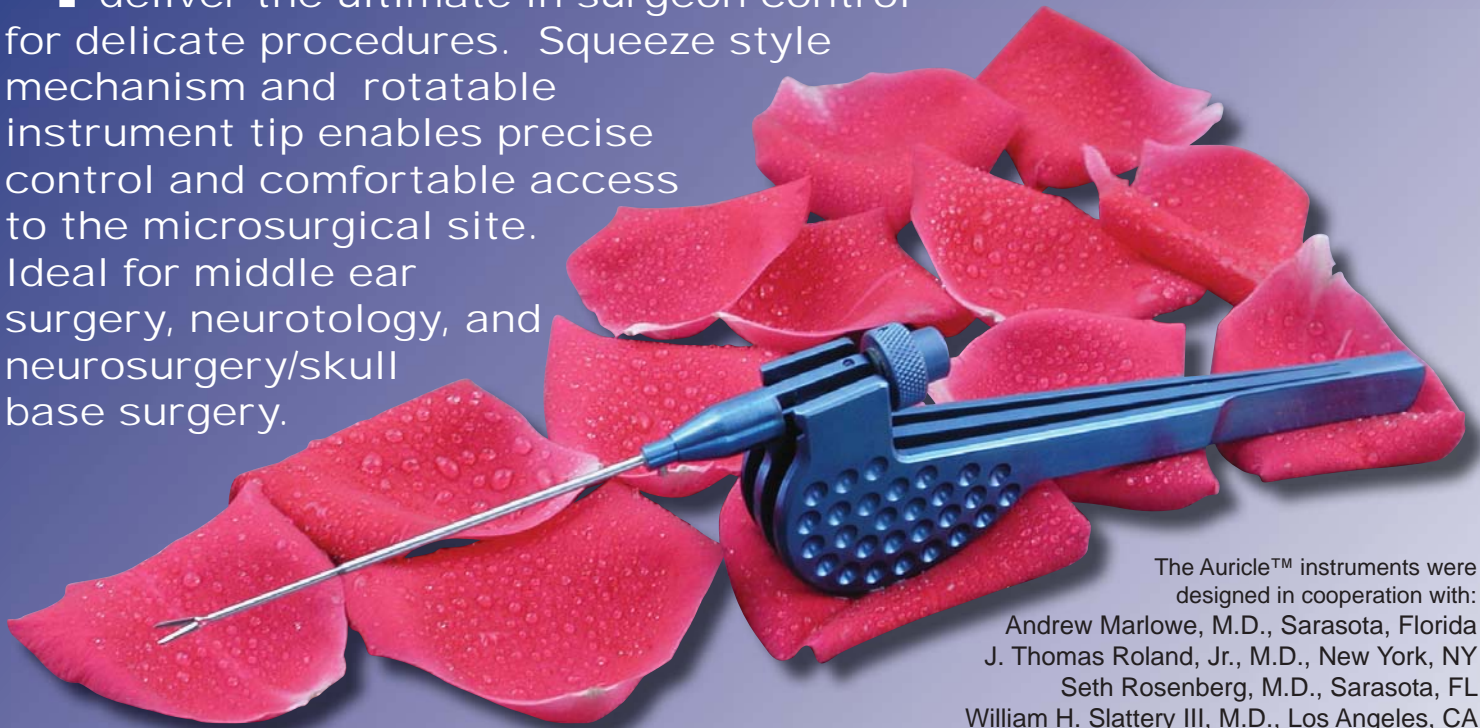

The Auricle™ instruments were designed in cooperation with:
Andrew Marlowe, M.D., Sarasota, Florida
J. Thomas Roland, Jr., M.D., New York, NY
Seth Rosenberg, M.D., Sarasota, FL
William H. Slattery III, M.D., Los Angeles, CA

- Unique squeeze style handle enables precise control of the instrument tip during actuation.
- Lightweight titanium handle with large gripping area for comfortable hand position reduces hand fatigue.
- Instrument tip rotates 360° allowing you to move the tip instead of your hand.
- Instrument shaft is .90mm allowing excellent visualization and easy access for working in small spaces.
- Shaft length of 55mm or 80mm available for various surgical approaches.
- A variety of unique scissors and forceps are available for specific applications.




MedOne
www.MedOne.com

Auricle™ Instruments

VISIT OUR WEBSITE TO SEE A VIDEO
OF THE AURICLE INSTRUMENTS



SCISSORS

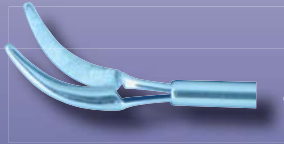

4001

4001 AURICLE CURVED SCISSORS (55MM)



4002

4002 AURICLE CURVED SCISSORS (80MM)

Blunt tipped curved blades.



4004

4002 AURICLE ANGLED SCISSORS 135° (55MM)



4005

4008 AURICLE ANGLED SCISSORS 135° (80MM)

Blunt tipped scissors angled 135°, 3mm from the tip.



4000

4004 AURICLE STRAIGHT SCISSORS (55MM)



4003

4010 AURICLE STRAIGHT SCISSORS (80MM)

Blunt straight scissors with 4mm long blades.



4040

4005 AURICLE SHARP TIP SCISSORS (55MM)



4042

4011 AURICLE SHARP TIP SCISSORS (80MM)

Pointed, straight blades 4mm long.

FORCEPS



4000

4000 AURICLE ANGLED FORCEPS (55MM)



4003

4006 AURICLE ANGLED FORCEPS (80MM)

Angled 2.5mm from the tip with slightly roughened jaws.



4040

4003 AURICLE SERRATED FORCEPS (55MM)



4042

4009 AURICLE SERRATED FORCEPS (80MM)

Blunt alligator style forceps with 4mm long jaws for general grasping. Jaws have serrations for non-slip grip.

PROBES



4040

4040 BALL END PROBE 90°*



4042

4041 BALL END PROBE 45°*

Thin .90mm shaft 80mm long, angled 3mm from tip.

SPATULA



4042

4042 SPATULA DISSECTOR*

Thin .90mm shaft 80mm long, angled 55°, 4mm from the tip.

*DEVELOPED IN COOPERATION WITH
J. THOMAS ROLAND, JR., M.D.,
NEW YORK, NEW YORK